



SAFE life™

T-Swab Cotine (COT) & Alcohol Oral Fluid Test

What does SAFElife™ T-Swab Cotine (COT) & Alcohol Oral Fluid Test do?

SAFElife™ T-Swab Cotine (COT) & Alcohol Oral Fluid Test is a rapid oral fluid test for screening exposure to tobacco smoke including cigarettes and vapor as well as alcohol using.

COT & Alcohol Detection

Cotine (COT), a biomarker for tobacco exposure detection, is an alkaloid found in tobacco and is also a major metabolite of Nicotine, which produces stimulation of the autonomic ganglia and central nervous system in humans.

Saliva alcohol test reflect the Blood Alcohol Contraction (BAC), the United States Department of Transportation (DOT) has established a BAC of 0.02% (0.02g/dL) as the cut off level at which an individual is considered positive for the presence of alcohol.

Optional Alcohol test

Ordering our product with or without Alcohol test upon your demand.
Optional Alcohol test with transparent window design for clear viewing.

Choose SAFElife™ T-Swab Cotine (COT) & Alcohol Oral Fluid Test for the detection of recent Nicotine use and Alcohol use.

- Light weighted pen-shaped device
- No privacy room needed for sample collection
- Simple, rapid, minimally-invasive saliva collection




Workflow

Ready to use & Soon to get results

Remove the cap



Place the collection pad under the tongue

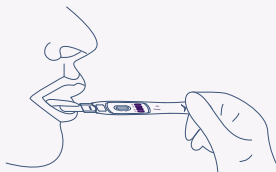

Get reliable results at 5 minutes for COT test,
at 2 minutes for Alcohol test



5 minutes



* Please refer to the Package Insert for detailed procedure.

Item code	Product name	Description	Complexity	Product image
SODOA-116-COT-P	SAFElife™ T-Swab Cotine (COT) Oral Fluid Test	1-Panel T-Swab; COT 30 ng/ml	Prescription Use Only	
SODOA-126NT-A	SAFElife™ T-Swab Cotine (COT) & Alcohol Oral Fluid Test	2-Panel T-Swab; COT+ALC	Forensic Use Only	

Rel: 2024/05/17